



ControMed

LED operating otoscope


Examination and treatment with one otoscope

- Compact & light weight.
- Wide opening in the back of the earpiece makes it easy to insert a surgical instrument.
- Earpiece is easily detached with one click for washing and sterilizing.



LED operating otoscope

Novel design enables both examination and treatment.

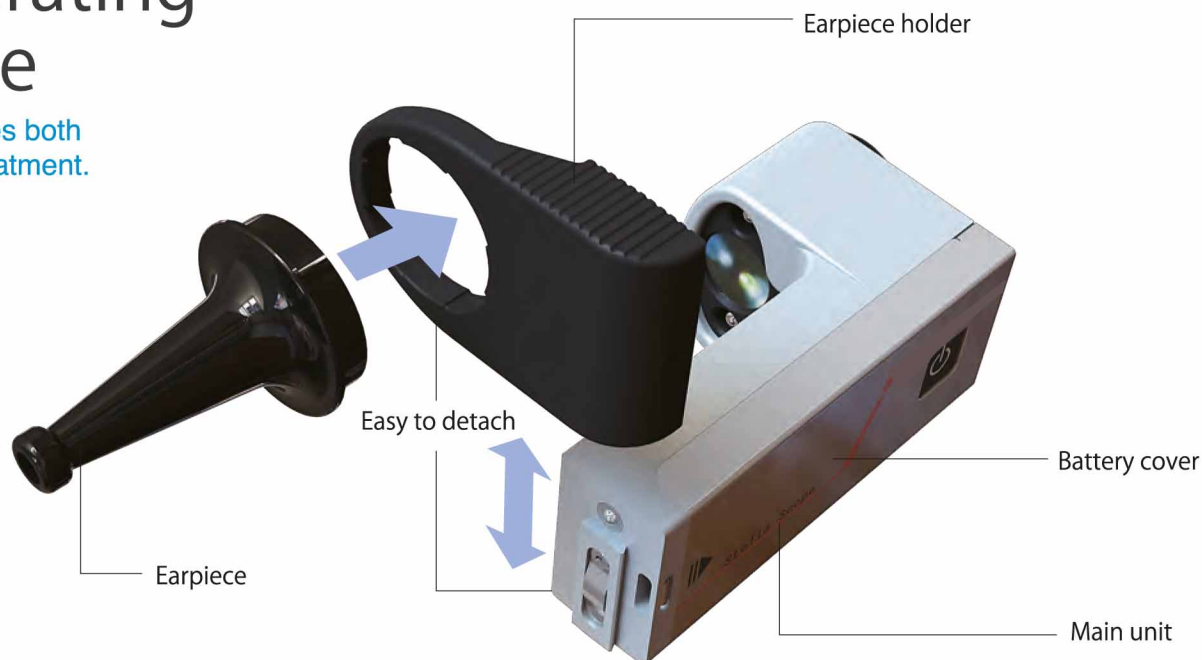

Features



Warm white LEDs with special optical design

Warm white light LEDs reveal the ear in familiar colors. This leads to safer treatment. Special optical design lets light penetrate deep into the ear.



Light weight and compact

The slim light weight design is comfortable to hold for a long time. The robust aluminum construction will last for a long time.



Wide opening in the back of the earpiece

Wide opening in the back of the earpiece makes it easy to insert a surgical instrument.

Easy maintenance earpiece and earpiece holder

The earpiece and the earpiece holder can be detached easily for washing and sterilizing.

8-hour battery operation

Two AAA batteries provide up to 8 hours operation. Rechargeable batteries can also be used.

Specifications

Dimensions	48 (W) x 97 (D) x 36 (H) mm
Weight	127 g
Power source	AAA battery x 2
Operating time	Approx. 8 hours
Light source	3 white LEDs

Composition

Main unit	1		
Earpiece holder (90-020-10)	1		
AAA battery	2		
Earpiece No. 00 (90-031-02)	1	Internal tip diameter	2.8 mm
Earpiece No. 0 (90-031-01)	1	~	3.2 mm
Earpiece No. 1 (90-031-10)	1	~	4.0 mm
Earpiece No. 2 (90-031-20)	1	~	5.0 mm
Earpiece No. 3 (90-031-30)	Option	~	6.8 mm
Earpiece No. 4 (90-031-40)	Option	~	8.0 mm

ControMed

POX-BOX: 1056, Doncaster Heights,
VIC 3109, Australia
Tel : +61 4 3074 2079
Fax: +61 3 9891 7050
Website: www.contromed.com